




Skin System

Vascular
Blemish
treatment
system

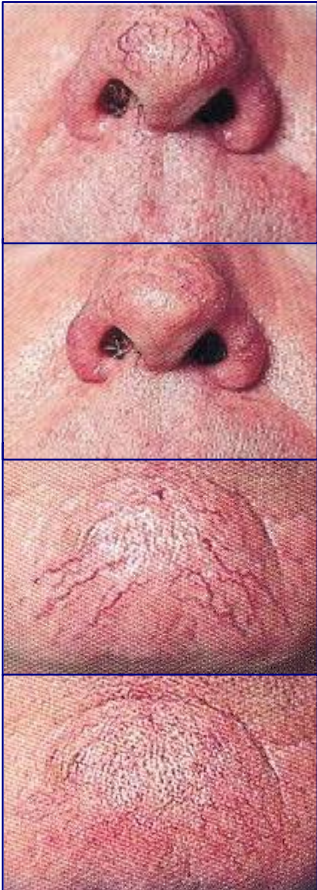

For quick and easy treatment of telangiectasia, cherry angiomas, Skin Tags, and Vascular blemishes.

One or two foot pedal operation for total control of each current

Low and High Frequency for telangiectasia and Skin tags

Phoresis – post treatment (cataphoresis) – Reduces redness, swelling, promotes healing

Air – for desensitizing the skin

4 pre-set keys for Small, Medium and Large capillaries

1 custom key to store a favorite setting

Test key for testing accessories for continuity



Ryan Masa – Representative

480-664-9490 or Toll Free 1-866-410-5549

www.skintechsystems.com

ryan@skintechsystems.com



Skin System – Vascular blemish treatment system

The Skin System is a quick and effective treatment for telangiectasia and skin tags.

The Skin System provides an inexpensive procedure for the treatment of Telangiectasia, cherry angiomas, spider veins, and skin tags found on the face and body. The Skin System provides better results than high intensity devices such as laser or hyfrecators.

1. Little discomfort
2. No bruising
3. More affordable



The quick and simple treatment will not cause scarring or hyperpigmentation. The Skin System uses low level Radio Frequency (RF) and DC (galvanic) current simultaneously to coagulate the capillary without desiccating surrounding tissue. The treatment is vascular specific resulting in minimal skin reaction.

The Phoresis modality will reduce redness; swelling and the germicidal effect of the current promotes the healing process.

The Skin System includes:

- 1 – Treatment cord
- 2 – Foot switches
- 1 – Stainless steel ground electrode
- 1 – Phoresis sphere electrode
- 1 – Set of accessory cords
- 1 – Power cord
- 1 – Procedure book “Telangiectasia” by Bono
- 1 – Skin System User Manual
- 50 – Treatment probes
- 3 year warranty

100-120 volts AC, 50/60 Hz

How long does a treatment take?

Most blemishes can be treated in a few seconds. Depending on the amount of telangiectasia grouping the treatment session can last about 10-15 minutes.

Is it painful?

Some experience a prickly sensation that can be mildly uncomfortable. Your skin care therapist may use a numbing agent for the treatment in sensitive areas to minimize the discomfort. If an injectable is used the formulation will be without epinephrine in order to diminish the effect of vessel constriction.

How is the procedure done?

Your skin care therapist will clean and prep the area. A sterile probe is inserted into the capillary and RF and DC current pass through the probe into the capillary and coagulates the blood flow. The coagulated segment is the clot that blocks the flow of blood in the capillary. When telangiectases are removed, new normal-size capillaries regrow to efficiently feed the skin.

Post treatment is critical for success. Nothing must be done to open the coagulated blood vessel segments. Anything that increases the skin’s blood flow or stretches the skin must be avoided for at least 24 hours. Clients must not stretch their skin, sunbathe, exercise, take a sauna,, wash their face or apply make-up for at least 24 hours following treatment.

Optional Upgrade!

Complement the *Skin System*.

Add

Red

and

Blue



Light Therapy to promote healthy skin